

**DIMAR**

PRODUCT  
LINE

Ambulance  
Emergency

# *DimAir Line*




***EASY VENT EMERGENCY  
MASK***

# EASY VENT EMERGENCY MASK

**CPAP "E.V.E. Mask" system for the administration of air/oxygen mixtures in high flow CPAP techniques on hypoxemic patients with cardiogenic edema or acute respiratory failure (ARF). Equipped with venturi device capable of delivering up to 120 Lt/min. Light and practical kit, usable in Ambulance, First Aid and in hospital wards, in presence of a single oxygen cylinder or flowmeter. Flow regulation determines the value of PEEP administered to the patient.**

## TECHNICAL FEATURES:

- High efficiency integrated single-use venturi.
- Single patient.
- Mask with adjustable anatomical cushion.
- Integrate manometer for constant control and regulation of the applied pressure.
- Luer M fitting for oxygen delivery or EtCO<sub>2</sub> control
- Oxygen supply tube of 2 meters
- Adjustable headgear, quick application and fixing
- Latex and DEHP Free

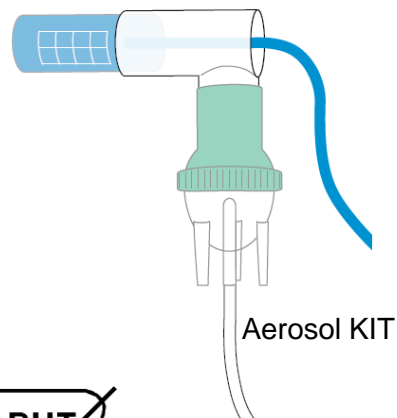

Aerosol KIT



## ADVANTAGES:

- Extremely Fast, Practical and Effective
- Continuous monitoring of applied pressure even in the presence of leaks
- The system's effectiveness allows for pressure up to approx. 20 cmH<sub>2</sub>O
- Reduced oxygen consumption for longer life of the cylinder
- For its functioning only one oxygen source is sufficient, like flow meter or cylinder
- Possibility of administering aerosols by special kit (code 700/5066)

## REFERENCES:

CODE	SIZE	RDM	CND/EMDN
700/6413	PED/ADULT (S/M)	2683358	R0301010202
700/6414	ADULT (L/XL)		

- In single pack.
- Possibility to create models, versions and custom packaging and on specific request.

- Multiple pack of 10 pieces.

## PERFORMANCE TABLE

Flow/Pressures/FiO<sub>2</sub> Flow obtained 60 Lt

Flow Lt./min	PEEP approx	O <sub>2</sub> Con. approx
5	5	40%
6	7	40%
7	9	40%
8	11	40%
9	14	40%
10	16	40%
11	19	40%



**HMC Premedical S.p.A**

Via Bosco 1/3 – 41037 Mirandola – MO – I  
Tel: +39 0535 22704 Fax: +39 0535 609546

